
100x140mm

Heated Massage Support Brace USER MANUAL

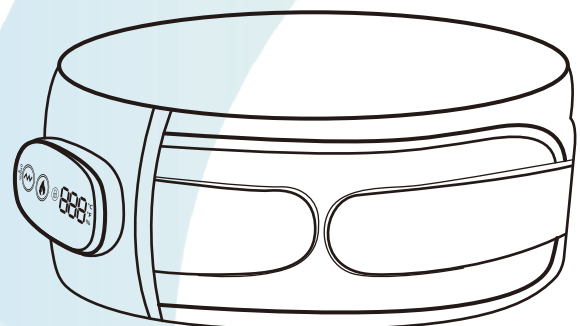MODEL: 102319

Thank you for choosing our product.
Please read this manual carefully before use. After reading,
please keep the manual properly for future reference.

CONTENT

EN	1-6
FR	7-12

PRODUCT DESCRIPTION

Welcome to choose our newly upgraded Heated Massage Support Brace. Our product boasts excellent elasticity and flexibility, enabling it to closely conform to the curves of your waist while providing personalized support and massage effects. Whether you are working, driving or doing other daily activities, it can effectively reduce waist pressure, improve waist problems, and provide you with a full range of waist care, so that you can enjoy intimate health care in your busy life.

WARM REMINDER: This product is for relief only and is not a substitute for medical equipment or medications. If any discomfort or abnormality occurs during use, please stop using it immediately and consult a professional.

PRODUCT PARAMETERS

Material	Neoprene
Battery Type	955565 polymer lithium battery
Battery Capacity	5000mAh/18.5wh
Working Voltage	5V
Charging Current	2000mA
Charging Time	3.5hours

EASY TO WEAR


USING STEPS

- 01** Please inspect the product before use to ensure there is no breakage or obvious damage. Choose a quiet, comfortable environment so you can completely relax and enjoy your massage.
- 02** Wear the product on the area to be used and adjust the position to ensure it fits tightly on the massage area.
- 03** Long press the power button of the controller to turn on the device. If it cannot be turned on normally, it may be due to insufficient battery. Please charge it and try again.
- 04** Touch the vibration button to turn on the massage function. The display area displays "1", indicating that the current vibration level is 1. Touch the vibration button repeatedly to switch vibration gears: 1-2-3-0. The display area shows "0" and the massage function is turned off.
* When using it for the first time, it is recommended to start with a lower gear and gradually increase the intensity after getting used to it.
- 05** Touch the heating button to turn on the heating mode, the display area shows "45 °C", and the current heating temperature level is 1. Touch the heating button repeatedly to switch the heating temperature: 45 °C-50 °C-55 °C-60 °C-65 °C-0. The display area shows "0" and the heating function is turned off.
* When using it for the first time, it is recommended to start with a lower temperature, gradually adapt to it, and then increase the temperature to avoid high temperature burns on the skin.
- 06** The controller is equipped with dual temperature units of °C and °F. Long press the heating button and vibration button at the same time to switch the temperature unit.

- 07** Short press the power button, the button lock light will light up red, the controller buttons are locked, and the display area will briefly display the current power. You cannot operate the controller at this time to avoid accidentally touching the buttons while using the product. To operate the controller, press the power button to unlock it.
- 08** This product has a 30-minute timer function. After 30 minutes of use, the device will automatically shut down.
- 09** During use, long press the power button until the display area displays "OFF" to turn off the device.
- 10** After use, store the product and accessories safely in a dry and ventilated place away from direct sunlight or moisture.

Note: If you feel uncomfortable or find any abnormalities in the product during use, please stop using it immediately and contact customer service personnel for processing.

CONTROLLER

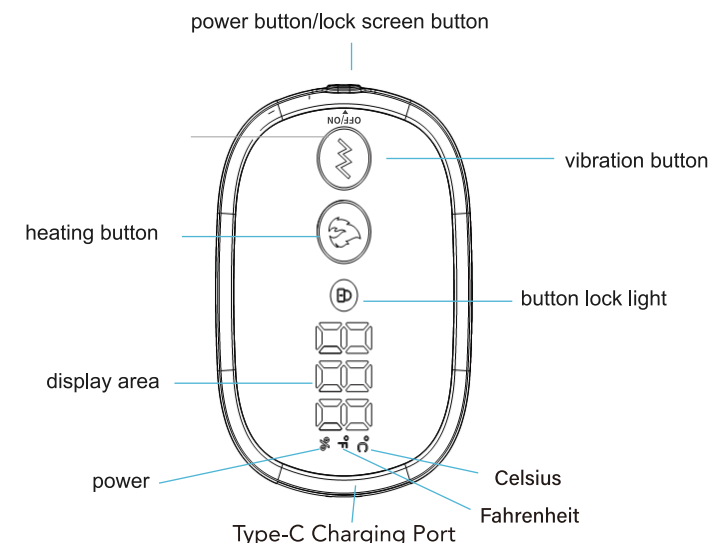Note: The numbers on the controller screen indicate the massage gear, heating temperature (°C/°F) and charging current power & progress (%).
1. When setting the heating function, the number + °C/°F represents the heating temperature.
2. When setting the vibration function, the number represents the massage gear.
3. When the product is in use (not charging), the number + % represents the current power.
4. When the product is charging, the number + % represents the charging progress.

ABOUT CHARGING

- When the device is powered on, short press the power button, and the display area will show the current battery level.
- When the screen shows that the power is less than 20%, please charge it in time.
- Please make sure the device is turned off before charging.
- Take out the equipped charging cable, connect to the controller and the charging adapter. The display screen of the controller lights up, indicating that the device starts charging, and the numbers represent the current power. During the charging process, the number will slowly increase.
- The controller display shows "100%", indicating that the device is fully charged.

Note: When charging, please use our original charging cable, otherwise there may be problems with of not charging or not charging fully. Do not operate the controller during charging, otherwise may damage it.

PRECAUTIONS


- Excessive massage is harmful to the body.**
 - This product is a household appliance and cannot replace medical equipment or medicines.**
 - Incorrect operation may result in safety accidents or property damage.**
- Please make sure the product is completely dry before use.
 - When starting to use this product, it is recommended to select the low temperature mode. Slight redness and warmth of the skin after use is normal. Please feel free to use it. However, if severe skin redness and swelling occur after use, please stop use immediately and consult a doctor.
 - The device is not waterproof and cannot be used in a wet place.
 - It is recommended to avoid direct contact of the product with the skin.
 - It is recommended to start with the lowest temperature and gradually adjust to the appropriate temperature according to personal conditions to avoid scalding the skin of people who are not sensitive to heat.
 - Long-term use will cause the product temperature to be high, so be careful of burns. (Continuous use of the heating function should not exceed 30 minutes)
 - If the user experiences any discomfort during use, such as dizziness, headache, asthma, etc., please stop using it immediately; if necessary, please contact your doctor.
 - When not in use, please turn off the product and keep it in a clean and dry place.
 - Do not expose the device to sunlight or fire to avoid damaging the thermal conductor of the insulator.
 - Do not disassemble or modify the controller and circuit board at will to avoid artificial damage.

WARNING

CONSULT A DOCTOR BEFORE USE:
● Patients suffering from osteoporosis, fractures, sprains, acute pain and other leg diseases
● Patients with nerve disorders in the legs or sensory disturbances in the skin
● Pregnant or postpartum women
● Heart disease patients
* The product has a heating function, people who are not sensitive to temperature must pay attention.

CLEANING AND MAINTENANCE

- Before cleaning, please make sure to disconnect the power completely.
- When cleaning, wipe the product with a damp cloth, but be careful not to wet the Type-C interface. Do not use strong cleaners or hard brushes.
- Lay flat to dry.
- This product does not support machine wash or hand wash.
- Do not dry clean.

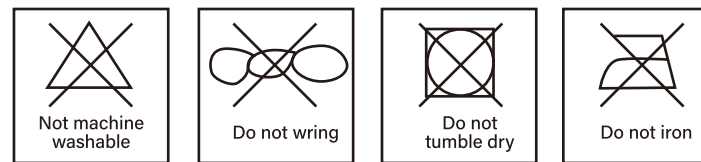WARRANTY & AFTER-SALES SERVICE

365 Days Quality Guarantee
Thank you for purchasing our product. Your satisfaction and support are our greatest motivation. Please contact us directly if you have any questions. We will provide prompt service if you encounter any problem with our product.
Thank you!

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